



WASHING AND DISINFECTION

Laboratory, pharma and industry sector

MORE THAN **70 YEARS** OF MADE IN ITALY



Founded in 1948 by the longest standing Italian family in the production of cooking appliances, Smeg is now the only appliance company with **PRODUCTION EXCLUSIVELY IN ITALY.**

The company name, in Italian, is an acronym of: Smalterie Metallurgiche Emiliane Guastalla (in English it translates as "Metal Enamelling Factory based in Guastalla, Emilia").



Products Categories



MEDICAL

Thermo-disinfectors and accessories for the medical field.



LABORATORY

Glassware washers for laboratories with a large choice in accessories to meet all washing requirements.



DENTAL

Comprehensive range of machines designed for the **automated cleaning in dental practices.**

As experts in glass washers and washers disinfectors, we control products and processes, we understand the market's varying needs and comply with the strictest international standards and directives.

Continuously we are monitoring and seeking to optimize our processes, products and performance.

Smeg GW Series: Under-Counter models

COMPLIANCE WITH THE FOLLOWING DIRECTIVES:

2006/42/EC

2014/30/UE

2011/65/UE

2015/863

COMPLIANCE WITH THE FOLLOWING STANDARDS:

EN 61010-1:2010

EN 61010-2-040:2015


EN 61326-1:2013

EN12100:2010

EN 50581:2012



LABORATORY

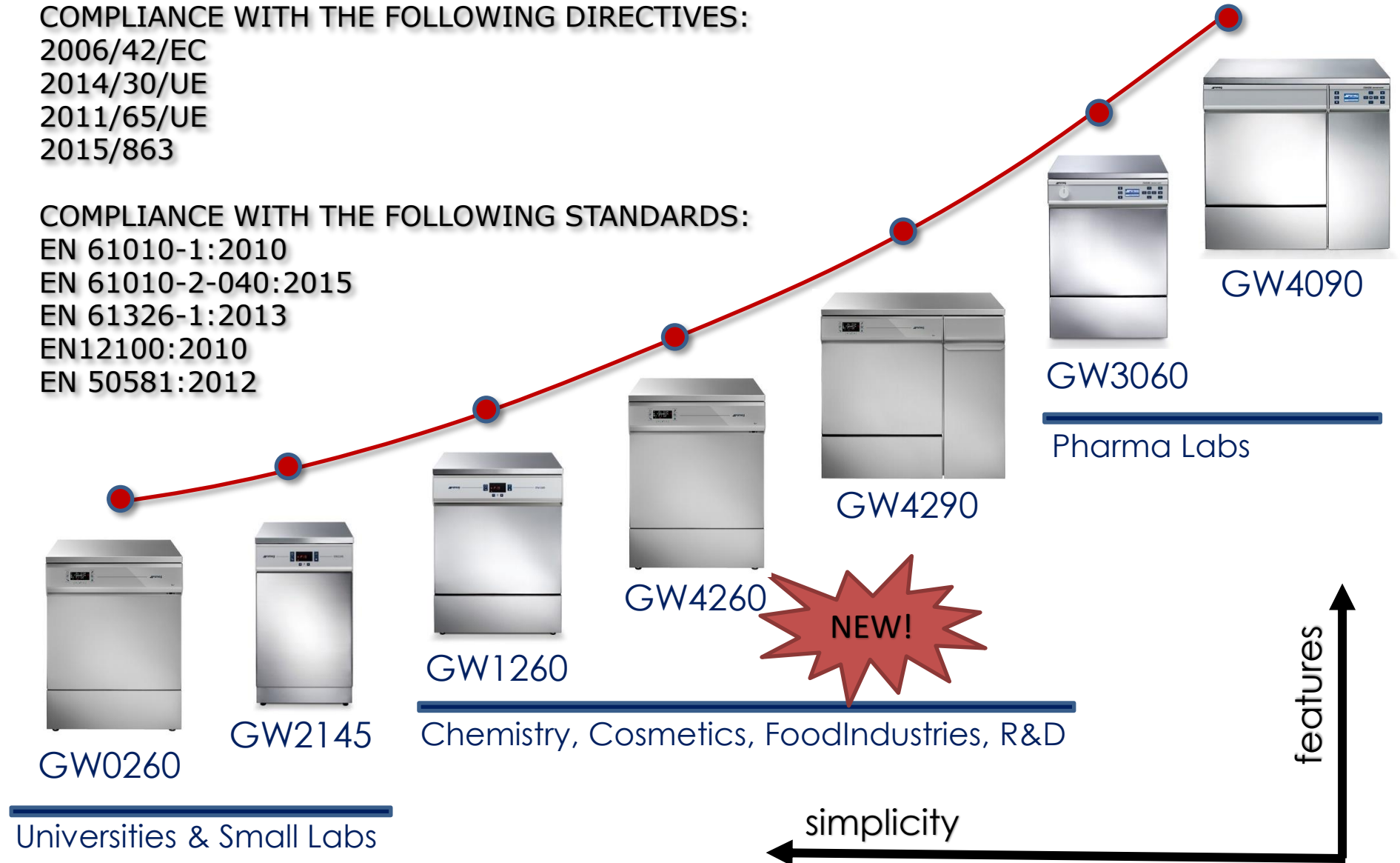

GlasswareWasher GW4260

The **GW4260** series is the result of more than 30 years of experience in washing and disinfection for laboratories. It combines the most advanced technologies and reliable solutions in full compliance with current directives and standards.

The innovative system constituted by glassware-washer and flexible accessories allow the reprocessing of a wide range of instruments as well as optimizing spaces and decreasing the end user's costs.

The electrical consumption has been widely reduced by combining the smart management of the electrical heating, a mindful design of washing cycles and an efficient active thermodynamic drying system.



LABORATORY



quality

sustainability

security

saving

Main Features

20 default Programs

20 custom Programs

Drying system included,
made of an hot air
generator (99.99% HEPA
filter)

software to set all the
parameters trough PC trough
USB connection

Real-time A0-value

Softener

Greater efficiency – treats
20% more liters of water

Glass Door available
as an option



All the undercounter models make
available a usable **washing area** up to
1m², depending on the rack.

The high precision in dosing of
detergents and flow meter control on
water intake minimize wastes as well as
the environmental impact is remarkably
reduced. Both the washing chamber
and the inner door are made of
stainless steel **AISI 316L** while the
external panels are made of stainless
steel AISI 304

Cooling water discharge
with mixed water, 65°C

Water inlet pipes with
Aquastop System

NEW!



LABORATORY

New Display




USB Connection
Included

RS232 port and LAN connection
are still available as options



LABORATORY

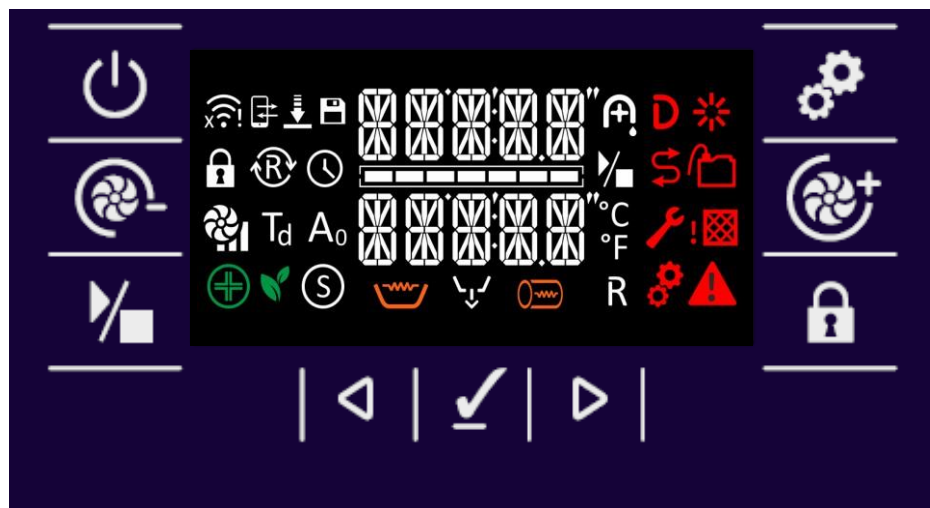

- Increased information available
- Program progress bar
- Intuitive icons
- Animated Icons
- Colored Icons

GW4260

New Codes




LABORATORY

VERSION	DETERGENT POWDER DISPENSER	NEUTRALIZER DOSING PUMP	DETERGENT DOSING PUMP	ADDITIVE DOSING PUMP	STEAM CONDENSER
GW4260-00-000	✓	✓	Not available	Optional	Not available
GW4260-0C-000	✓	✓	Not available	Optional	✓
GW4260-S0-000	Not available	✓	✓	Optional	Not available
GW4260-SC-000	Not available	✓	✓	Optional	✓

Dimensions



LABORATORY

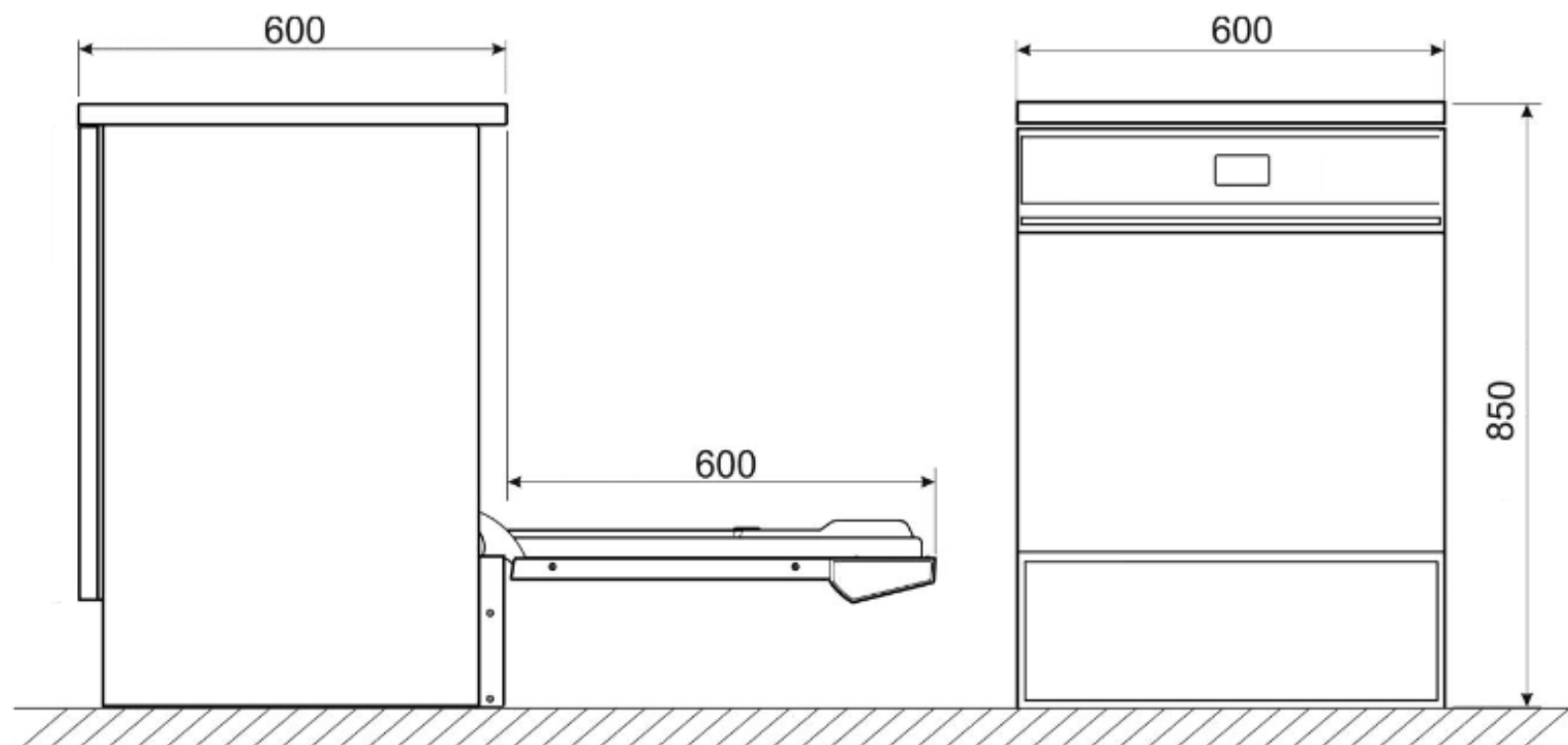

Accessories

TROLLEY FOR BUTYROMETERS
WITH DRYING SYSTEM



TROLLEY FOR
FISCHER TUBES



2 LEVELS INJECTION TROLLEY FOR BOTTLES



WASHING TROLLEY FOR FLASKS
PIPETTES / TEST TUBES - 20
POSITIONS



LABORATORY



www.smeg-instruments.com